

	Premarket Notification 510(k) - FDA	FORM. 02.1 - V	
		Issuance:	19/12/2023
		Revision:	00
		Sheet:	1/12

Instructions for Use

Stitch - Suture Wires and Tapes.

DEVICE DESCRIPTION

Stitch - Suture Wires and Tapes are used in surgeries for approximations and for union of soft tissues until healing occurs. The surgical sutures and tapes are made of UHMWPE (ultra-high molecular weight polyethylene), and the needles are made in AISI300 steel series, according to the standard ASTM F899. They are found in a variety of diameters to meet the range of anatomies of the patients; it is provided in sterile condition sterilized by Ethylene Oxide. The sutures are classified as non-absorbable and meet the USP standards. This product is for single use only. This device comes in sizes USP 3-0, 2-0, 0, 1, 2 and 5 and suture tape sizes 1.3mm, 1.6mm and 2.0mm.

The color additives used are accepted by FDA, according to the guidance "Summary of Color Additives for Use in the United States in Foods, Drugs, Cosmetics, and Medical Devices"; Color Additives Approved for Use in Medical Devices, as follows:

- Chromium-cobalt-aluminum oxide (CFR §73.1015) for the blue color in concentration lower than 2% by weight of the suture material.
- D&C Black (CFR §74.3054) for the black color in concentration lower than 1% by weight of the suture material.

INDICATIONS FOR USE

Stitch - Suture Wires and Tapes are indicated for use in general soft tissue approximation and/or ligation.

PERFORMANCE

The mechanical analysis shows that the tensile strength of the suture wires and tapes did not vary significantly during the shelf life of the product and all tests showed compliance with the mechanical specifications of the US Pharmacopeia. Also, the chemical analysis confirms the declared chemical composition and shows that the material did not degrade over time.

CONTRAINDICATIONS

Contraindications are listed below, but are not limited to:

- Patients with active infection, fragile tissue, or poor vascularity;
- Patients with illnesses or any clinical condition that may compromise the healing process.

Note: Patients who smoke, are elderly, and suffer from degenerative diseases or metabolic dysfunction (e.g., diabetes) have greater difficulty achieving tissue regeneration. Consequently, there is a greater chance of tissue regeneration/healing not occurring, which may lead to thread and/or tape breakage, separation of unhealed soft tissue, or infection.

	Premarket Notification 510(k) - FDA	FORM. 02.1 - V	
		Issuance:	19/12/2023
		Revision:	00
		Sheet:	2/12

WARNINGS AND PRECAUTIONS

- This device is intended to be used by a trained medical professional.
- Do not re-sterilize this device.
- Once open, discard unused suture.
- The sutures must never be re-used.
- This is a single use device. Reuse of this device could result in failure of the device to perform as intended and could cause harm to the patient and/or user.
- Detailed instructions on the use and limitations of this device should be given to the patient.
- In handling this or any other suture material, care should be taken to avoid damage from handling. Avoid crushing or crimping damage due to application of surgical instruments such as forceps or needle holders.
- Surgeons must apply their professional judgment when determining the appropriate sized device based on the specific indication, preferred surgical technique, and patient history.

POTENTIAL ADVERSE EVENTS

- Infections, both deep and superficial.
- Foreign body reactions.

HOW SUPPLIED

The Stitch family is supplied in pre-cut lengths, in various configurations, with needles, to meet the range of anatomies of the patients; it is provided in sterile condition sterilized by Ethylene Oxide.

MRI SAFETY INFORMATION

The implantable devices are made only in polymers with its main component the UHMWPE (Polyethylene Ultra-High Molecular Weight). The Stitch - Suture Wires and Tapes is MR Safe.

STORAGE

Sterile devices must be stored in the original unopened packaging, away from moisture and should not be used after the expiration date.

The hospital is responsible for the steps of handling, conservation, and storage of the products from receipt, and must ensure that the manufacturer's recommendations are respected, and that the product maintains the established original characteristics.

The product must be preserved, handled, and transported in such a way as to prevent any damage or alteration to its characteristics and packaging.

	Premarket Notification 510(k) - FDA	FORM. 02.1 - V	
		Issuance:	19/12/2023
		Revision:	00
		Sheet:	3/12

The product must be conserved and stored in its original packaging until the moment of use, with the due care of the medical-hospital area, in a clean, dry, airy environment, protected from sunlight and free of contaminating substances (acid and organic vapors), and with controls recommended, at room temperature. The handling of the product must be done exclusively by professionals in the medical-hospital area, duly trained, qualified and familiar with the technique and procedures involved.

SURGICAL TECHNIQUE

1. Carefully open the packaging and remove the product;
2. Apply the product in the desired region, in order to carry out the planned technique;
3. Repeat these steps as many times as necessary;
4. Cut the needle, when provided with the product, and any excess thread/ribbon, and discard;
5. Clean the place;
6. Close surgical access.

USER WARNING

These Instructions for Use are available on-line, through the electronic address of the manufacturer: <http://www.gmreis.com.br/produtos/IFU>, and can be checked in the search field by trade name, described on the label of the product packaging. The Instructions for use provided will always be in accordance with the current latest version. If there is interest from the user, the Instructions for Use may be available in printed format, without additional cost. Ask for free by e-mail: sac@gmreis.com.br

Caution: Federal (USA) law restricts this device to sale by or on the order of a licensed physician.

FURTHER INFORMATION

If the medical product presents a risk specific unpredictable, being outside of its specifications or being generated any dissatisfaction, directly notify the GM Dos Reis Customer Service.

GM Dos Reis Indústria e Comércio Ltda - Pierre Simon de Laplace Avenue, no.600 TECHNOPARK – CEP: 13069-320 – Campinas – SP – BRAZIL or directly notify at Telephone Number:+55 19 3765-9900; E-mail: sac@gmreis.com.br

SYMBOLS GLOSSARY

ANSI/AAMI/ISO 15223-1 and ASTM F2503-13.

Symbol	Description
	Manufacturing date
	Manufactured by
	Expiration date
	Keep dry
Rx only	Prescription use only
STERILE EO	Sterile by Ethylene Oxide
	Do not reuse
LOT	Lot number
	Caution "see instructions for use"
	Keep away from sunlight
	Do not use it if package is damaged
	Do not re-sterilize
COD	Product catalogue number
MR	MR Safe
QT	Quantity
UDI	Unique Device Identifier

COMPANY NAME / MANUFACTURER NAME:

GM dos Reis Industria e Comercio Ltda.

Avenida Pierre Simon de Laplace, no. 600 Lote 3 - Quadra F - Block 9677 TECHNOPARK – CEP: 13069-320 – Campinas – SP - BRAZIL

Email: qualidade@gmreis.com.br

Telephone number: +55 19 37659900

STITCH EN - Revision: 0, October 15, 2025

Premarket Notification 510(k) - FDA

FORM. 02.1 - V


Issuance: 19/12/2023

Revision: 00

Sheet: 5/12

The Stitch - Suture Wires and Tapes are available in the dimensions listed below.

STITCH

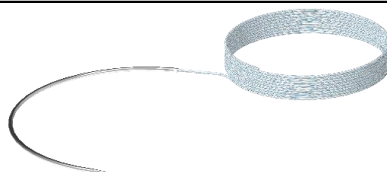

CODE	DESCRIPTION
333-30AB-3813CI	STITCH - UHMWPE Suture wire USP 3-0 45cm white and blue with needle 3/8 13mm cylindrical tip
333-30PB-3813CI	STITCH - UHMWPE Suture wire USP 3-0 45cm white and black with needle 3/8 13mm cylindrical tip
333-30B-3813CI	STITCH - UHMWPE Suture wire USP 3-0 45cm white with needle 3/8 13mm cylindrical tip
333-30A-3813CI	STITCH - UHMWPE Suture wire USP 3-0 45cm blue with needle 3/8 13mm cylindrical tip
333-30AB-3815CI	STITCH - UHMWPE Suture wire USP 3-0 45cm white and blue with needle 3/8 15mm cylindrical tip
333-30PB-3815CI	STITCH - UHMWPE Suture wire USP 3-0 45cm white and black with needle 3/8 15mm cylindrical tip
333-30B-3815CI	STITCH - UHMWPE Suture wire USP 3-0 45cm white with needle 3/8 15mm cylindrical tip
333-30A-3815CI	STITCH - UHMWPE Suture wire USP 3-0 45cm blue with needle 3/8 15mm cylindrical tip
333-30AB-3816R	STITCH - UHMWPE Suture wire USP 3-0 45cm white and blue with needle 3/8 16mm reverse tip
333-30PB-3816R	STITCH - UHMWPE Suture wire USP 3-0 45cm white and blue with needle 3/8 16mm reverse tip
333-30B-3816R	STITCH - UHMWPE Suture wire USP 3-0 45cm white with needle 3/8 16mm reverse tip
333-30A-3816R	STITCH - UHMWPE Suture wire USP 3-0 45cm blue with needle 3/8 16mm reverse tip
333-30AB-3826DI	STITCH - UHMWPE Suture wire USP 3-0 45cm white and blue with needle 3/8 26mm diamond tip
333-30PB-3826DI	STITCH - UHMWPE Suture wire USP 3-0 45cm white and black with needle 3/8 26mm diamond tip
333-30B-3826DI	STITCH - UHMWPE Suture wire USP 3-0 45cm white with needle 3/8 26mm diamond tip
333-30A-3826DI	STITCH - UHMWPE Suture wire USP 3-0 45cm blue with needle 3/8 26mm diamond tip
333-30AB-1213CI	STITCH - UHMWPE Suture wire USP 3-0 45cm white and blue with needle 1/2 13mm cylindrical tip
333-30PB-1213CI	STITCH - UHMWPE Suture wire USP 3-0 45cm white and black with needle 1/2 13mm cylindrical tip
333-30B-1213CI	STITCH - UHMWPE Suture wire USP 3-0 45cm white with needle 1/2 13mm cylindrical tip
333-30A-1213CI	STITCH - UHMWPE Suture wire USP 3-0 45cm blue with needle 1/2 13mm cylindrical tip
333-2-036-AB3817CI	STITCH - UHMWPE Suture wire USP 2-0 45cm blue and white with needle 3/8 17mm cylindrical tip
333-2-036-PB3817CI	STITCH - UHMWPE Suture wire USP 2-0 45cm white and black with needle 3/8 17mm cylindrical tip
333-2-036-B3817CI	STITCH - UHMWPE Suture wire USP 2-0 45cm white with needle 3/8 17mm cylindrical tip
333-2-036-A3817CI	STITCH - UHMWPE Suture wire USP 2-0 45cm blue with needle 3/8 17mm cylindrical tip

Premarket Notification 510(k) - FDA

FORM. 02.1 - V

Issuance: 19/12/2023

Revision: 00

Sheet: 6/12

333-2-090-AB3817CI	STITCH - UHMWPE Suture wire USP 2-0 90cm white and blue with needle 3/8 17mm cylindrical tip
333-2-090-PB3817CI	STITCH - UHMWPE Suture wire USP 2-0 90cm white and black with needle 3/8 17mm cylindrical tip
333-2-090-B3817CI	STITCH - UHMWPE Suture wire USP 2-0 90cm white with needle 3/8 17mm cylindrical tip
333-2-090-A3817CI	STITCH - UHMWPE Suture wire USP 2-0 90cm blue with needle 3/8 17mm cylindrical tip
333-2-018-AB1226CI	STITCH - UHMWPE Suture wire USP 2-0 45cm white and blue with needle 1/2 26mm cylindrical tip
333-2-018-PB1226CI	STITCH - UHMWPE Suture wire USP 2-0 45cm white and black with needle 1/2 26mm cylindrical tip
333-2-018-B1226CI	STITCH - UHMWPE Suture wire USP 2-0 45cm white with needle 1/2 26mm cylindrical tip
333-2-018-A1226CI	STITCH - UHMWPE Suture wire USP 2-0 45cm blue with needle 1/2 26mm cylindrical tip
333-2-090-AB1226CI	STITCH - UHMWPE Suture wire USP 2-0 90cm white and blue with needle 1/2 26mm cylindrical tip
333-2-090-PB1226CI	STITCH - UHMWPE Suture wire USP 2-0 90cm white and black with needle 1/2 26mm cylindrical tip
333-2-090-B1226CI	STITCH - UHMWPE Suture wire USP 2-0 90cm white with needle 1/2 26mm cylindrical tip
333-2-090-A1226CI	STITCH - UHMWPE Suture wire USP 2-0 90cm blue with needle 1/2 26mm cylindrical tip
333-0-AB-1226CI	STITCH - UHMWPE Suture wire USP 0 90cm white and blue with needle 1/2 26mm cylindrical tip
333-0-PB-1226CI	STITCH - UHMWPE Suture wire USP 0 90cm white and black with needle 1/2 26mm cylindrical tip
333-0-B-1226CI	STITCH - UHMWPE Suture wire USP 0 90cm white with needle 1/2 26mm cylindrical tip
333-0-A-1226CI	STITCH - UHMWPE Suture wire USP 0 90cm blue with needle 1/2 26mm cylindrical tip
333-0-AB-1226DI	STITCH - UHMWPE Suture wire USP 0 90cm white and blue with needle 1/2 26mm diamond tip
333-0-PB-1226DI	STITCH - UHMWPE Suture wire USP 0 90cm white and black with needle 1/2 26mm diamond tip
333-0-B-1226DI	STITCH - UHMWPE Suture wire USP 0 90cm white with needle 1/2 26mm diamond tip
333-0-A-1226DI	STITCH - UHMWPE Suture wire USP 0 90cm blue with needle 1/2 26mm diamond tip
333-190-AB12265CI	STITCH - UHMWPE Suture wire USP 1 90cm white and blue with needle 1/2 26mm cylindrical tip
333-190-PB12265CI	STITCH - UHMWPE Suture wire USP 1 90cm white and black with needle 1/2 26mm cylindrical tip
333-190-B1226CI	STITCH - UHMWPE Suture wire USP 1 90cm white with needle 1/2 26mm cylindrical tip
333-190-A1226CI	STITCH - UHMWPE Suture wire USP 1 90cm blue with needle 1/2 26mm cylindrical tip
333-239-AB1226CI	STITCH - UHMWPE Suture wire USP 2 90cm blue and white with needle 1/2 26mm cylindrical tip
333-239-PB1226CI	STITCH - UHMWPE Suture wire USP 2 90cm white and black with needle 1/2 26mm cylindrical tip
333-239-B1226CI	STITCH - UHMWPE Suture wire USP 2 90cm white with needle 1/2 26mm cylindrical tip
333-239-A1226CI	STITCH - UHMWPE Suture wire USP 2 90cm blue with needle 1/2 26mm cylindrical tip
333-290-AB12265R	STITCH - UHMWPE Suture wire USP 2 90cm white and blue with needle 1/2 26mm reverse tip
333-290-PB12265R	STITCH - UHMWPE Suture wire USP 2 90cm white and black with needle 1/2 26mm reverse tip
333-290-B12265R	STITCH - UHMWPE Suture wire USP 2 90cm white with needle 1/2 26mm reverse tip
333-290-A12265R	STITCH - UHMWPE Suture wire USP 2 90cm blue with needle 1/2 26mm reverse tip

Premarket Notification 510(k) - FDA


FORM. 02.1 - V

Issuance: 19/12/2023

Revision: 00

Sheet: 7/12

333-290-AB12366CI	STITCH - UHMWPE Suture wire USP 2 90cm white and blue with needle ½ 36mm cylindrical tip
333-290-PB12366CI	STITCH - UHMWPE Suture wire USP 2 90cm white and black with needle ½ 36mm cylindrical tip
333-290-B1236CI	STITCH - UHMWPE Suture wire USP 2 90cm white with needle ½ 36mm cylindrical tip
333-290-A1236CI	STITCH - UHMWPE Suture wire USP 2 90cm blue with needle ½ 36mm cylindrical tip
333-290-AB12366R	STITCH - UHMWPE Suture wire USP 2 90cm white and blue with needle ½ 36mm reverse tip
333-290-PB12366R	STITCH - UHMWPE Suture wire USP 2 90cm white and black with needle ½ 36mm reverse tip
333-290-B1236R	STITCH - UHMWPE Suture wire USP 2 90cm white with needle ½ 36mm reverse tip
333-290-A1236R	STITCH - UHMWPE Suture wire USP 2 90cm blue with needle ½ 36mm reverse tip
333-536-AB1248CO	STITCH – UHMWPE Suture wire USP 5 90cm white and blue with needle ½ 48mm conventional tip
333-536-PB1248CO	STITCH – UHMWPE Suture wire USP 5 90cm white and black with needle ½ 48mm conventional tip
333-536-B1248CO	STITCH – UHMWPE Suture wire USP 5 90cm white with needle ½ 48mm conventional tip
333-536-A1248CO	STITCH – UHMWPE Suture wire USP 5 90cm blue with needle ½ 48mm conventional tip
333-590-AB1248CI	STITCH – UHMWPE Suture wire USP 5 90cm white and blue with needle ½ 48mm cylindrical tip
333-590-PB1248CI	STITCH – UHMWPE Suture wire USP 5 90cm white and black with needle ½ 48mm cylindrical tip
333-590-B1248CI	STITCH – UHMWPE Suture wire USP 5 90cm white with needle ½ 48mm cylindrical tip
333-590-A1248CI	STITCH – UHMWPE Suture wire USP 5 90cm blue with needle ½ 48mm cylindrical tip

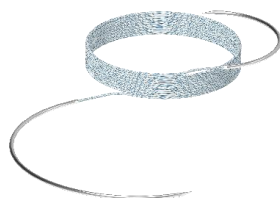

CODE	DESCRIPTION
333-30AB2-3813CI	STITCH - UHMWPE Suture wire USP 3-0 45cm white and blue with 2 needles 3/8 13mm cylindrical
333-30PB2-3813CI	STITCH - UHMWPE Suture wire USP 3-0 45cm white and black with 2 needles 3/8 13mm cylindrical
333-30B-23813CI	STITCH - UHMWPE Suture wire USP 3-0 45cm white with 2 needles 3/8 13mm cylindrical tip
333-30A-23813CI	STITCH - UHMWPE Suture wire USP 3-0 45cm blue with 2 needles 3/8 13mm cylindrical tip
333-30AB2-3815CI	STITCH - UHMWPE Suture wire USP 3-0 45cm white and blue with 2 needles 3/8 15mm cylindrical
333-30PB2-3815CI	STITCH - UHMWPE Suture wire USP 3-0 45cm white and black with 2 needles 3/8 15mm cylindrical
333-30B-23815CI	STITCH - UHMWPE Suture wire USP 3-0 45cm white with 2 needles 3/8 15mm cylindrical tip
333-30A-23815CI	STITCH - UHMWPE Suture wire USP 3-0 45cm blue with 2 needles 3/8 15mm cylindrical tip
333-30AB2-3816R	STITCH - UHMWPE Suture wire USP 3-0 45cm white and blue with 2 needles 3/8 16mm reverse

Premarket Notification 510(k) - FDA

FORM. 02.1 - V

Issuance: 19/12/2023

Revision: 00

Sheet: 8/12

333-30PB2-3816R	STITCH - UHMWPE Suture wire USP 3-0 45cm white and black with 2 needles 3/8 16mm reverse
333-30B-23816R	STITCH - UHMWPE Suture wire USP 3-0 45cm white with 2 needles 3/8 16mm reverse tip
333-30A-23816R	STITCH - UHMWPE Suture wire USP 3-0 45cm blue with 2 needles 3/8 16mm reverse tip
333-30AB2-3826DI	STITCH - UHMWPE Suture wire USP 3-0 45cm white and blue with 2 needles 3/8 26mm diamond tip
333-30PB2-3826DI	STITCH - UHMWPE Suture wire USP 3-0 45cm white and black with 2 needles 3/8 26mm diamond tip
333-30B-23826DI	STITCH - UHMWPE Suture wire USP 3-0 45cm white with 2 needles 3/8 26mm diamond tip
333-30A-23826DI	STITCH - UHMWPE Suture wire USP 3-0 45cm blue with 2 needles 3/8 26mm diamond tip
333-30AB2-1213CI	STITCH - UHMWPE Suture wire USP 3-0 45cm white and blue with 2 needles 1/2 13mm cylindrical tip
333-30PB2-1213CI	STITCH - UHMWPE Suture wire USP 3-0 45cm white and black with 2 needles 1/2 13mm cylindrical tip
333-30B-21213CI	STITCH - UHMWPE Suture wire USP 3-0 45cm white with 2 needles 1/2 13mm cylindrical tip
333-30A-21213CI	STITCH - UHMWPE Suture wire USP 3-0 45cm blue with 2 needles 1/2 13mm cylindrical tip
333-2045-2AB3817CI	STITCH - UHMWPE Suture wire USP 2-0 45cm white and blue with 2 needles 3/8 17mm cylindrical tip
333-2045-2PB3817CI	STITCH - UHMWPE Suture wire USP 2-0 45cm white and black with 2 needles 3/8 17mm cylindrical tip
333-2045-2B3817CI	STITCH - UHMWPE Suture wire USP 2-0 45cm white with 2 needles 3/8 17mm cylindrical tip
333-2045-2A3817CI	STITCH - UHMWPE Suture wire USP 2-0 45cm blue with 2 needles 3/8 17mm cylindrical tip
333-2090-2AB3817CI	STITCH - UHMWPE Suture wire USP 2-0 90cm white and blue with 2 needles 3/8 17mm cylindrical tip
333-2090-2PB3817CI	STITCH - UHMWPE Suture wire USP 2-0 90cm white and black with 2 needles 3/8 17mm cylindrical tip
333-2090-2B3817CI	STITCH - UHMWPE Suture wire USP 2-0 90cm white with 2 needles 3/8 17mm cylindrical tip
333-2090-2A3817CI	STITCH - UHMWPE Suture wire USP 2-0 90cm blue with 2 needles 3/8 17mm cylindrical tip
333-20-2AB1226CI	STITCH - UHMWPE Suture wire USP 2-0 45cm white and blue with 2 needles 1/2 26mm cylindrical tip
333-20-2PB1226CI	STITCH - UHMWPE Suture wire USP 2-0 45cm white and black with 2 needles 1/2 26mm cylindrical tip
333-20-2B1226CI	STITCH - UHMWPE Suture wire USP 2-0 45cm white with 2 needles 1/2 26mm cylindrical tip
333-20-2A1226CI	STITCH - UHMWPE Suture wire USP 2-0 45cm blue with 2 needles 1/2 26mm cylindrical tip
333-2090-2AB1226CI	STITCH - UHMWPE Suture wire USP 2-0 90cm white and blue with 2 needles 1/2 26mm cylindrical tip
333-2090-2PB1226CI	STITCH - UHMWPE Suture wire USP 2-0 90cm white and black with 2 needles 1/2 26mm cylindrical tip
333-2090-2B1226CI	STITCH - UHMWPE Suture wire USP 2-0 90cm white with 2 needles 1/2 26mm cylindrical tip
333-2090-2A1226CI	STITCH - UHMWPE Suture wire USP 2-0 90cm blue with 2 needles 1/2 26mm cylindrical tip
333-0-AB2-1226CI	STITCH - UHMWPE Suture wire USP 0 90cm white and blue with 2 needles 1/2 26mm cylindrical tip
333-0-PB2-1226CI	STITCH - UHMWPE Suture wire USP 0 90cm white and black with 2 needles 1/2 26mm cylindrical tip
333-0-B2-1226CI	STITCH - UHMWPE Suture wire USP 0 90cm white with 2 needles 1/2 26mm cylindrical tip
333-0-A2-1226CI	STITCH - UHMWPE Suture wire USP 0 90cm blue with 2 needles 1/2 26mm cylindrical tip
333-0-AB2-1226DI	STITCH - UHMWPE Suture wire USP 0 90cm white and blue with 2 needles 1/2 26mm diamond tip

Premarket Notification 510(k) - FDA

FORM. 02.1 - V

Issuance: 19/12/2023

Revision: 00

Sheet: 9/12

333-0-PB2-1226DI	STITCH - UHMWPE Suture wire USP 0 90cm white and black with 2 needles 1/2 26mm diamond tip
333-0-B2-1226DI	STITCH - UHMWPE Suture wire USP 0 90cm white with 2 needles 1/2 26mm diamond tip
333-0-A2-1226DI	STITCH - UHMWPE Suture wire USP 0 90cm blue with 2 needles 1/2 26mm diamond tip
333-190-2B1226CI	STITCH - UHMWPE Suture wire USP 1 90cm white with 2 needles 1/2 26mm cylindrical tip
333-190-2AB1226CI	STITCH - UHMWPE Suture wire USP 1 90cm white and blue with 2 needles 1/2 26mm cylindrical tip
333-190-2PB1226CI	STITCH - UHMWPE Suture wire USP 1 90cm white and black with 2 needles 1/2 26mm cylindrical tip
333-190-2A1226CI	STITCH - UHMWPE Suture wire USP 1 90cm blue with 2 needles 1/2 26mm cylindrical tip
333-290-AB12265CI	STITCH - UHMWPE Suture wire USP 2 90cm white and blue with 2 needles 1/2 26mm cylindrical tip
333-290-PB12265CI	STITCH - UHMWPE Suture wire USP 2 90cm white and black with 2 needles 1/2 26mm cylindrical tip
333-290-B21226CI	STITCH - UHMWPE Suture wire USP 2 90cm white with 2 needles 1/2 26mm cylindrical tip
333-290-A21226CI	STITCH - UHMWPE Suture wire USP 2 90cm blue with 2 needles 1/2 26mm cylindrical tip
333-290-2B12265R	STITCH - UHMWPE Suture wire USP 2 90cm white with 2 needles 1/2 26mm reverse tip
333-290-2B1236CI	STITCH - UHMWPE Suture wire USP 2 90cm white with 2 needles 1/2 36mm cylindrical tip
333-290-2B1236R	STITCH - UHMWPE Suture wire USP 2 90cm white with 2 needles 1/2 36mm reverse tip
333-536-2B1248CO	STITCH – UHMWPE Suture wire USP 5 90cm white with 2 needles 1/2 48mm conventional tip
333-590-2B1248CI	STITCH – UHMWPE Suture wire USP 5 90cm white with 2 needles 1/2 48mm cylindrical tip
333-590-AB1248CC	STITCH – UHMWPE Suture wire USP 5 90cm white and blue with needle 1/2 48mm conventional tip and needle 1/2 26mm cylindrical
333-590-BP1248CC	STITCH – UHMWPE Suture wire USP 5 90cm white and black with needle 1/2 48mm conventional tip and needle 1/2 26mm cylindrical
333-590-B1248CC	STITCH – UHMWPE Suture wire USP 5 90cm white with needle 1/2 48mm conventional tip and needle 1/2 26mm cylindrical
333-590-A1248CC	STITCH – UHMWPE Suture wire USP 5 90cm blue with needle 1/2 48mm conventional tip and needle 1/2 26mm cylindrical



CODE	DESCRIPTION
333-290-ABR76DI	STITCH - UHMWPE Suture wire USP 2 90cm white and blue with straight needle and diamond tip
333-290-PBR76DI	STITCH - UHMWPE Suture wire USP 2 90cm white and black with straight needle and diamond tip
333-290-BR76DI	STITCH - UHMWPE Suture wire USP 2 90cm white with straight needle and diamond tip

Premarket Notification 510(k) - FDA


FORM. 02.1 - V

Issuance: 19/12/2023

Revision: 00

Sheet: 10/12

333-290-AZ76DI	STITCH - UHMWPE Suture wire USP 2 90cm blue with straight needle and diamond tip
----------------	--

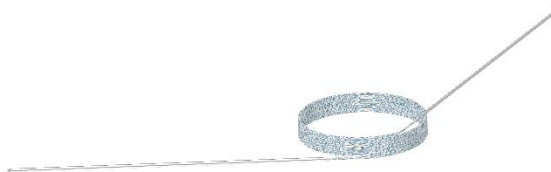


CODE	DESCRIPTION
333-290-AB2R76DI	STITCH - UHMWPE Suture wire USP 2 90cm white and blue with 2 straight needles and diamond tip
333-290-PB2R76DI	STITCH - UHMWPE Suture wire USP 2 90cm white and black with 2 straight needles and diamond tip
333-290-B2R76DI	STITCH - UHMWPE Suture wire USP 2 90cm white with 2 straight needles and diamond tip
333-290-A2R76DI	STITCH - UHMWPE Suture wire USP 2 90cm blue with 2 straight needles and diamond tip

STITCH LOOP STRAIGHT



CODE	DESCRIPTION
333-STL-2ABRTCI	STITCH LOOP STRAIGHT - UHMWPE Suture wire USP 2 with 50cm white and blue with straight needle 76,0mm cylindrical tip
333-STL-2PBRTCI	STITCH LOOP STRAIGHT - UHMWPE Suture wire USP 2 with 50cm white and black with straight needle 76,0mm cylindrical tip
333-STL-2BRTCI	STITCH LOOP STRAIGHT - UHMWPE Suture wire USP 2 with 50cm white with straight needle 76,0mm cylindrical tip
333-STL-2ARTCI	STITCH LOOP STRAIGHT - UHMWPE Suture wire USP 2 with 50cm blue with straight needle 76,0mm cylindrical tip

STITCH TAPE

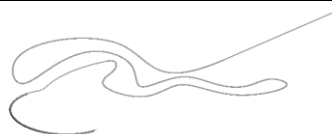

CODE	DESCRIPTION
333-13-B	STITCH TAPE - UHMWPE Tape 100cm white and black with 1,3mm with needle ½ 26,0mm cylindrical tip
333-13-AB	STITCH TAPE - UHMWPE Tape 100cm white and blue with 1,3mm with needle ½ 26,0mm cylindrical tip

Premarket Notification 510(k) - FDA

FORM. 02.1 - V

Issuance: 19/12/2023

Revision: 00

Sheet: 11/12

333-13-BR	STITCH TAPE - UHMWPE Tape 100cm white 1,3mm with needle ½ 26,0mm cylindrical tip
333-13-A	STITCH TAPE - UHMWPE Tape 100cm blue with 1,3mm with needle ½ 26,0mm cylindrical tip
333-16-B	STITCH TAPE - UHMWPE Tape 100cm white and black with 1,6mm with needle ½ 26,0mm cylindrical tip
333-16-AB	STITCH TAPE - UHMWPE Tape 100cm white and blue with 1,6mm with needle ½ 26,0mm cylindrical tip
333-16-BR	STITCH TAPE - UHMWPE Tape 100cm white 1,6mm with needle ½ 26,0mm cylindrical tip
333-16-A	STITCH TAPE - UHMWPE Tape 100cm blue with 1,6mm with needle ½ 26,0mm cylindrical tip
333-20-B	UHMWPE Tape 100cm white and black with 2,0mm with needle ½ 26,0mm cylindrical tip
333-20-AB	UHMWPE Tape 100cm blue and white 2,0mm with needle ½ 26,0mm cylindrical tip
333-20-BR	STITCH TAPE - UHMWPE Tape 100cm white 2,0mm with needle ½ 26,0mm cylindrical tip
333-20-A	STITCH TAPE - UHMWPE Tape 100cm blue with 2,0mm with needle ½ 26,0mm cylindrical tip
333-FTL13613-BP	STITCH LONG TAPE - UHMWPE Tape 136cm white and black with 1,3mm with needle ½ 26,0mm cylindrical tip
333-FTL13613-BA	STITCH LONG TAPE - UHMWPE Tape 136cm white and blue with 1,3mm with needle ½ 26,0mm cylindrical tip
333-FTL13613-BR	STITCH LONG TAPE - UHMWPE Tape 136cm white 1,3mm with needle ½ 26,0mm cylindrical tip
333-FTL13613-A	STITCH LONG TAPE - UHMWPE Tape 136cm blue with 1,3mm with needle ½ 26,0mm cylindrical tip
333-FTL13616-BP	STITCH LONG TAPE - UHMWPE Tape 136cm white and black with 1,6mm with needle ½ 26,0mm cylindrical tip
333-FTL13616-BA	STITCH LONG TAPE - UHMWPE Tape 136cm white and blue with 1,6mm with needle ½ 26,0mm cylindrical tip
333-FTL13616-BR	STITCH LONG TAPE - UHMWPE Tape 136cm white 1,6mm with needle ½ 26,0mm cylindrical tip
333-FTL13616-A	STITCH LONG TAPE - UHMWPE Tape 136cm blue with 1,6mm with needle ½ 26,0mm cylindrical tip
333-FTL1362-BP	STITCH LONG TAPE - UHMWPE Tape 136cm white and black with 2,0mm with needle ½ 26,0mm cylindrical tip
333-FTL1362-BA	STITCH LONG TAPE - UHMWPE Tape 136cm white and blue with 2,0mm with needle ½ 26,0mm cylindrical tip
333-FTL1362-BR	STITCH LONG TAPE - UHMWPE Tape 136cm white 2,0mm with needle ½ 26,0mm cylindrical tip
333-FTL1362-A	STITCH LONG TAPE - UHMWPE Tape 136cm blue with 2,0mm with needle ½ 26,0mm cylindrical tip

STITCH LACE



CODE	DESCRIPTION
333-ECAB-2-1248RB	STITCH LACE - UHMWPE Lace suture wire USP 2 with 66cm white and blue reinforced with closed loop 30.0mm with ½ 26.0mm needle, cylindrical tip.
333-ECB-2-1248RB	STITCH LACE - UHMWPE Lace suture wire USP 2 with 66cm white reinforced with closed loop 30.0mm with ½ 26.0mm needle, cylindrical tip.
333-ECBP-2-1248RB	STITCH LACE - UHMWPE Lace suture wire USP 2 with 66cm white and black reinforced with closed loop 30.0mm with ½ 26.0mm needle, cylindrical tip.




Premarket Notification 510(k) - FDA

FORM. 02.1 - V

Issuance: 19/12/2023

Revision: 00

Sheet: 12/12

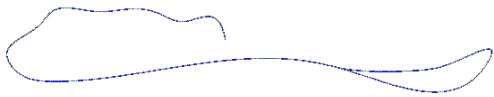333-ECA-2-1248RB	STITCH LACE - UHMWPE Lace suture wire USP 2 with 66cm blue reinforced with closed loop 30.0mm with ½ 26.0mm needle, cylindrical tip.
STITCH LINK	
	
CODE	DESCRIPTION
333-2AB-661226CI	STITCH LINK - UHMWPE Suture wire USP 2 with 66cm white and blue with closed loop with 135mm with ½ 26.0mm needle, cylindrical tip.
333-2B-661226CI	STITCH LINK - UHMWPE Suture wire USP 2 with 66cm white with closed loop with 135mm with ½ 26.0mm needle, cylindrical tip.
333-2BP-661226CI	STITCH LINK - UHMWPE Suture wire USP 2 with 66cm white and black with closed loop with 135mm with ½ 26.0mm needle, cylindrical tip.
333-2A-661226CI	STITCH LINK - UHMWPE Suture wire USP 2 with 66cm blue with closed loop with 135mm with ½ 26.0mm needle, cylindrical tip.